

Spinsolve 80




Experience the most powerful benchtop NMR spectrometer



Outstanding Features


- Highest field strength (80 MHz)
- Largest chemical shift spread
- Highest sensitivity
- Superb Resolution (0.5 Hz/20 Hz)
- Multi-nuclear probe
- 5 mm standard NMR tubes
- No sample spinning required
- Advanced methods like COSY, HSQC, HMBC, etc.
- Benchtop size and weight
- No cryogenics

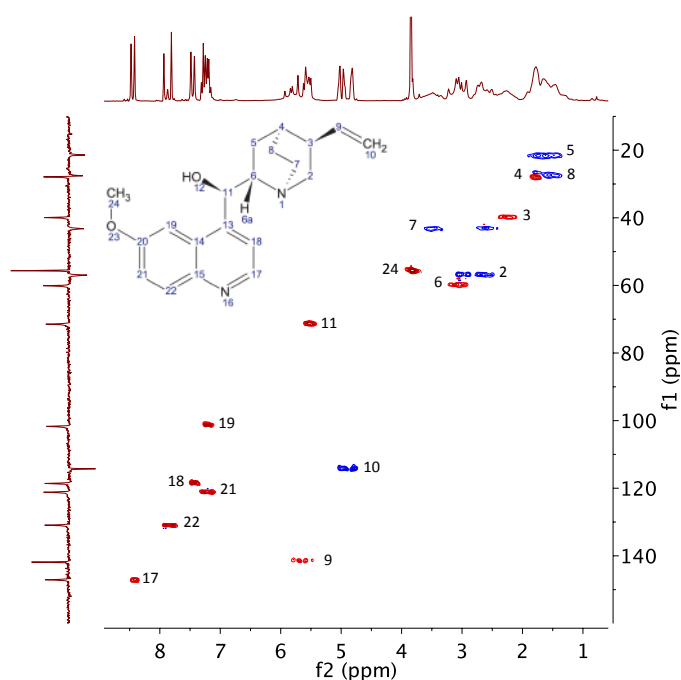


80 MHz HSQC-ME of Quinine

The most advanced multi-nuclear methods... in one click!


80 MHz NMR spectra of Ibuprofen

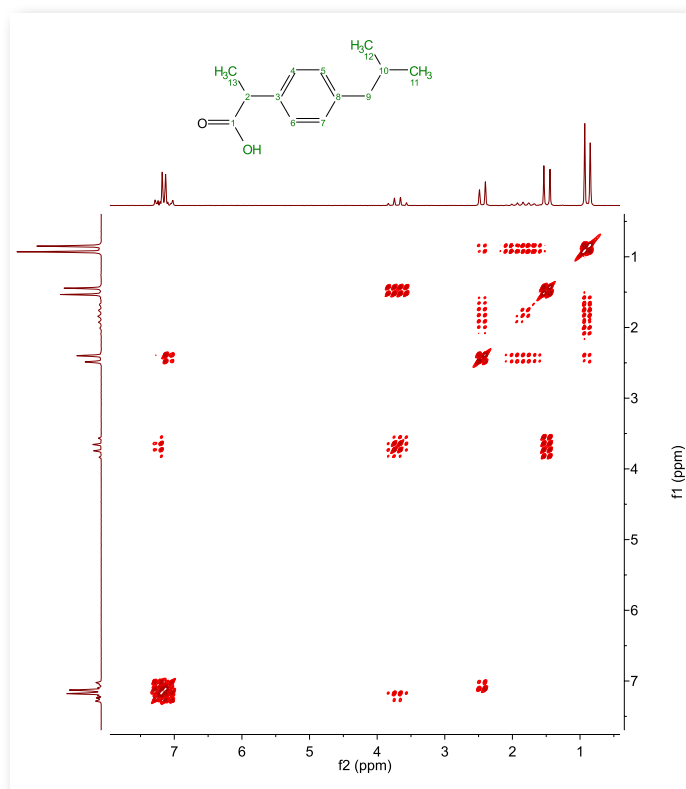


H-H COSY

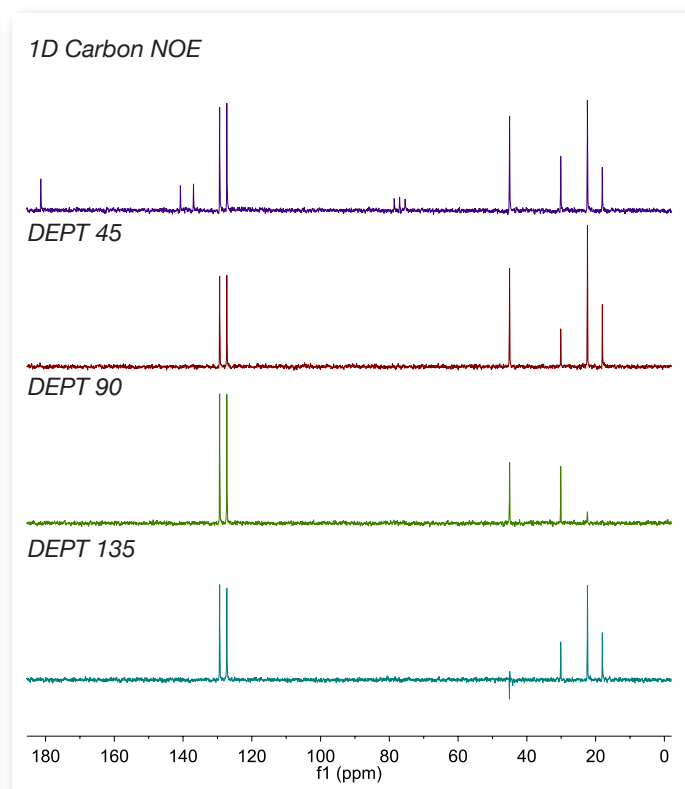

1D Carbon plus DEPT 45, 90, and 135 spectra

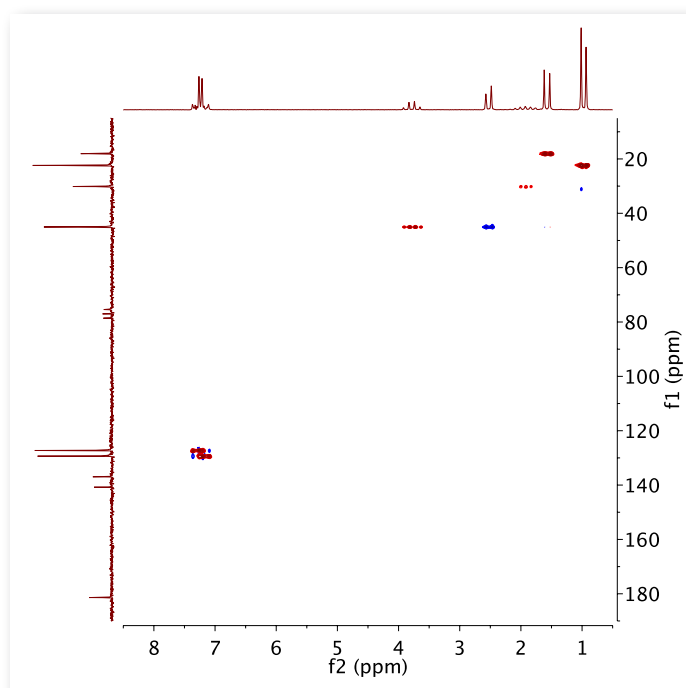

HSQC-ME

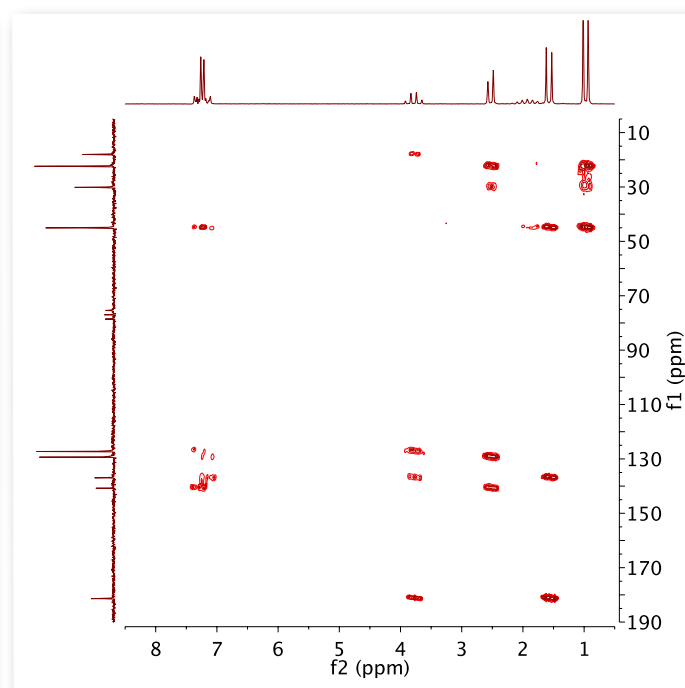

HMBC


^{13}C spectroscopy exploited for structure elucidation

Cholesterol

C₂₇H₄₆O, Mw: 386.65 g/mol

Concentration: 400 mMolar

Solvent: CDCl₃

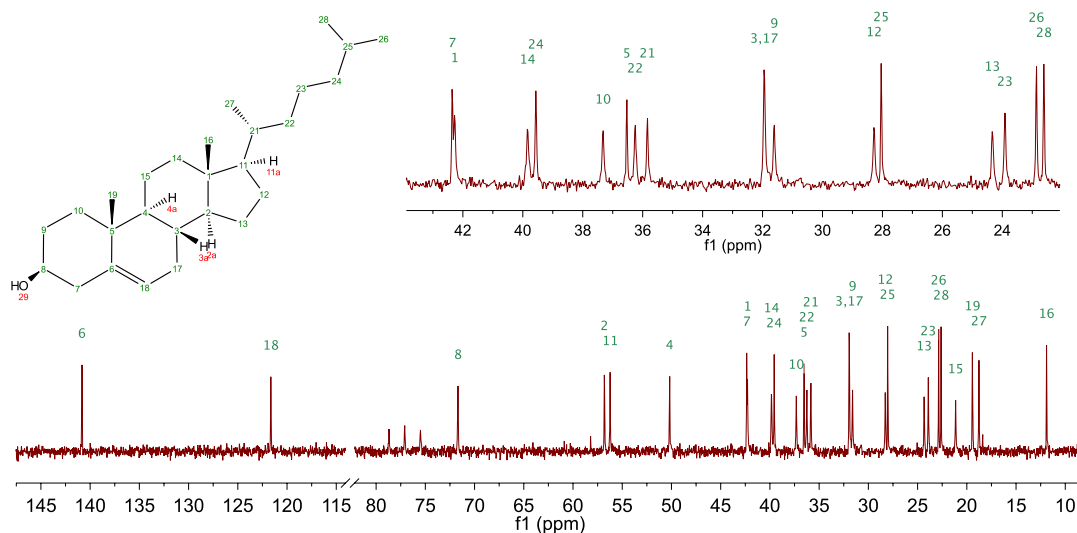


Hydro-Cortisone

C₂₁H₃₀O₅, Mw: 362.4 g/mol

Concentration: 500 mMolar

Solvent: DMSO

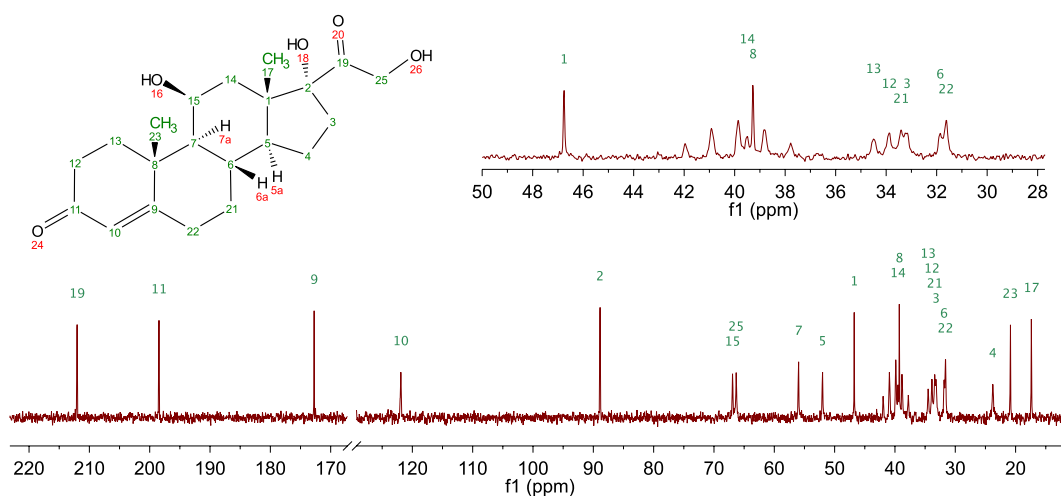


alpha-Linolenic

C₁₈H₃₀O₂, Mw: 278.44 g/mol

Concentration: 1 Molar

Solvent: CDCl₃

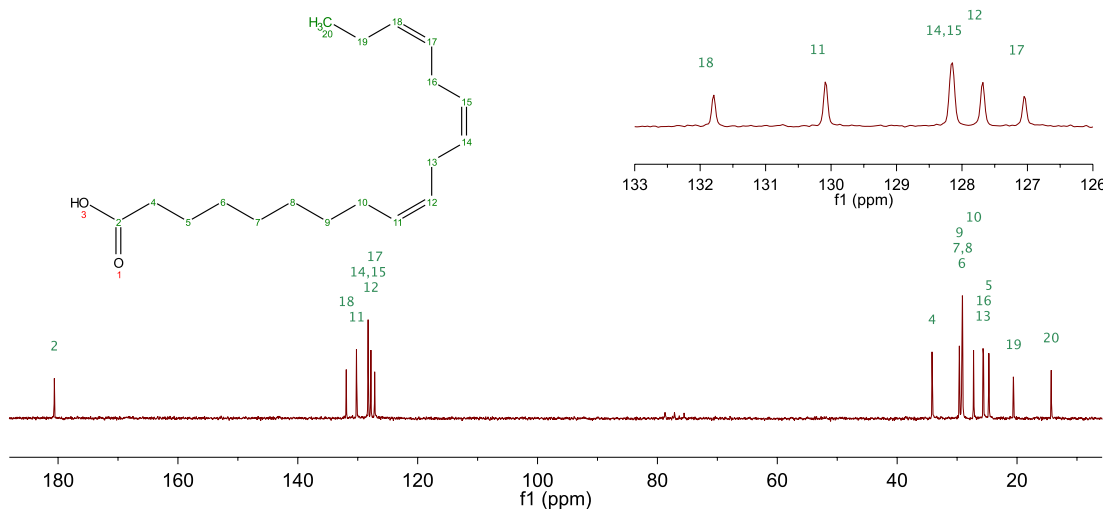



Specifications

- Nuclei: ^1H , ^{19}F , ^{13}C (Other Nuclei Available)
- Operating frequency: 80 MHz (^1H)
- ^1H 50% Linewidth: < 0.5 Hz
- ^1H 0.55% Linewidth: < 20 Hz
- ^1H Sensitivity (dual channel): >200:1 for 1% Ethyl Benzene
- Operating Temperature Range: 20° C to 25° C (68° F to 77° F)
- Dimensions: 58 x 43 x 40 cm (23" x 17" x 16")
- Weight: 72.5 kg (160 lb)
- Stray Field: < 2 G all around system
- Voltage Requirement: 100-240 VAC, 50/60 Hz



Pulse sequences available on the Spinsolve 80 Carbon spectrometer

Proton	Fluorine	Carbon
1D	1D	1D
Paramagnetic	Paramagnetic	DEPT
2D COSY	2D F - COSY	HETCOR
2D TOCSY	2D F - JRES	HMBC
2D JRES	2D FH - COSY	HMQC
T_1 , T_2	Reaction Monitoring	HSQC
Reaction Monitoring		HSQC-ME

Other sequences available, contact Magritek for details.

Contact us now for a quote, to request a demo or to measure your samples

CONTACT INFORMATION

For further information, please contact: sales@magritek.com

GERMANY

Philipsstraße 8
52068 Aachen, Germany
Tel: +49 (241) 70525-6000
Fax: +49 (241) 539-4375

NEW ZEALAND

6 Hurring Place, Unit 3
Newlands, Wellington 6037, NZ
Tel: +64 4 477 7096
Fax: +64 4 471 4665

UNITED STATES

6440 Lusk Blvd (D108)
San Diego, CA 92121, USA
Tel: +1 (855) 667-6835
+1 (866) NMR-MTEK

Or visit our website www.magritek.com